

Specifications

Power Supply	100~240VAC, 50/60Hz, 200VA
Flow Rate	0.1~10ml/s
Pressure Limit	325PSI
Auto Filling Rate	8ml/s
Air Purge Rate	A: 1.0~9.9ml/s B: 2.5~9.9ml/s
Piston Rod Speed	9.9ml/s
Injection Delay	1~300s
Scan Delay	1~300s
Pause Phase	1~900s
Phases	6
Injection Programs	100
Injection Records	100
Volume	A: 1~65ml
	B: 1~115ml
Multiple Language User Interface	



DC Power

High Quality Disposables



65ml Syringe for Contrast
115ml Syringe for Saline

antmed

Shenzhen Antmed Co., Ltd.

18 Jinhui Ave., Pingshan New District, Shenzhen 518122, China Tel: +86 755 8606 0992 Fax: +86 755 8606 0991 E-mail: info@antmed.com http://www.antmed.com





MR Injector

Leading Contrast Media Injectors Forward



MR Injector

Introducing an all-new ImaStar MR Injector, offering an extended list of product features, reliable performance, and a safe patient experience.




Ease of Operation and Clinical Efficiency

- · Fill contrast and saline with a single press of a button
- Optional IV pole

Independent Dual Motor Drive Systems

- Dual motor drive capability for time saving loading and filling of syringes
- · Allows independent injection of contrast and saline

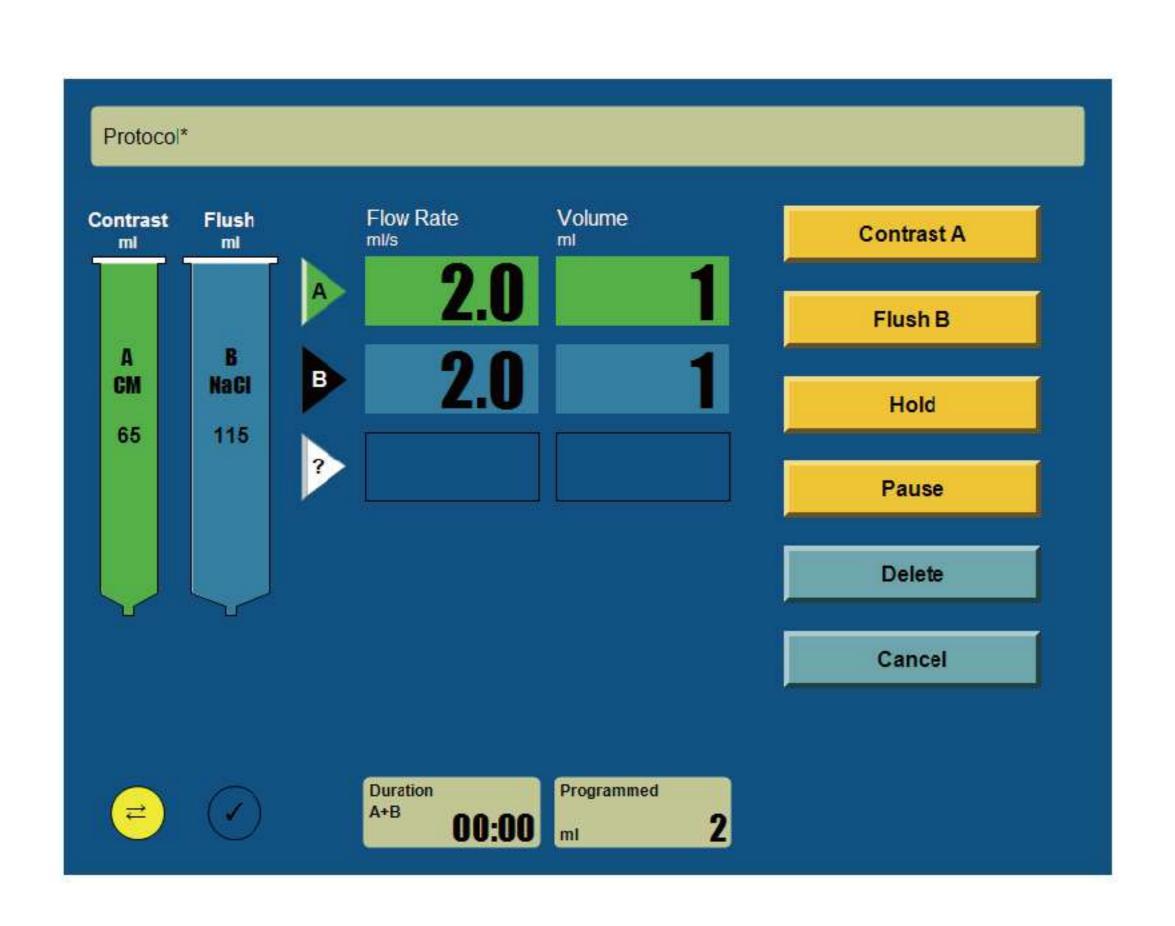



Color Touch-Screen

- · 12.1"color touch-screens for easy operation
- Clear and concise information display to ensure user operational confidence


Injection Record Storage

· Stores up to 100 injection history records




Real-time Pressure Curve

 Display provides real-time pressure graph for injection monitoring

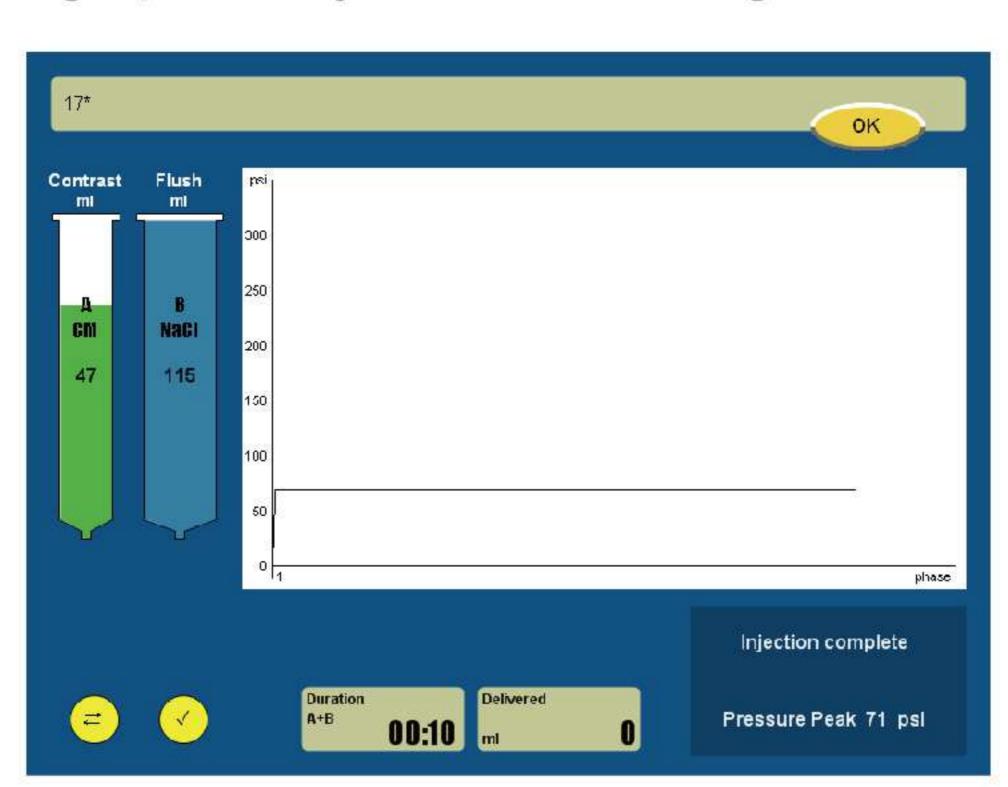


A Variety of Configuration Options

- Stores up to 100 customized injection programs
- Program storage capacity can be increased based on user needs

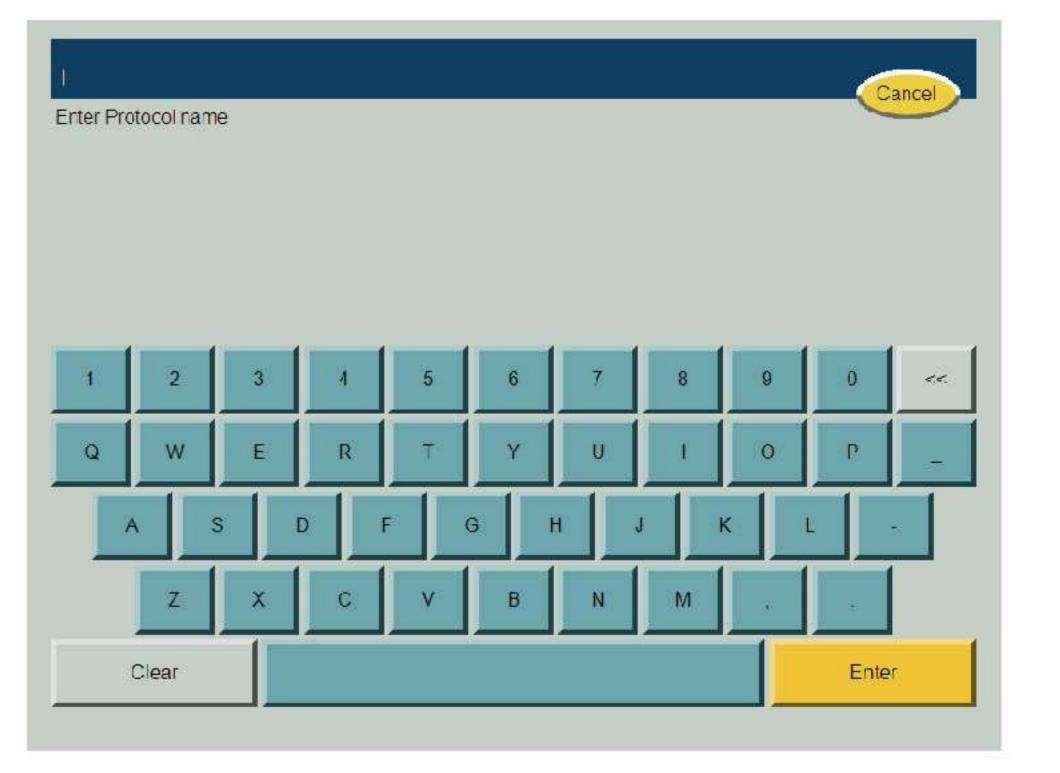

KVO Function

- Independent KVO (Keep Vein Open) ensures efficient examination procedure
- Flexible KVO interval time settings
- 115ml saline syringe capacity allows for longer KVO

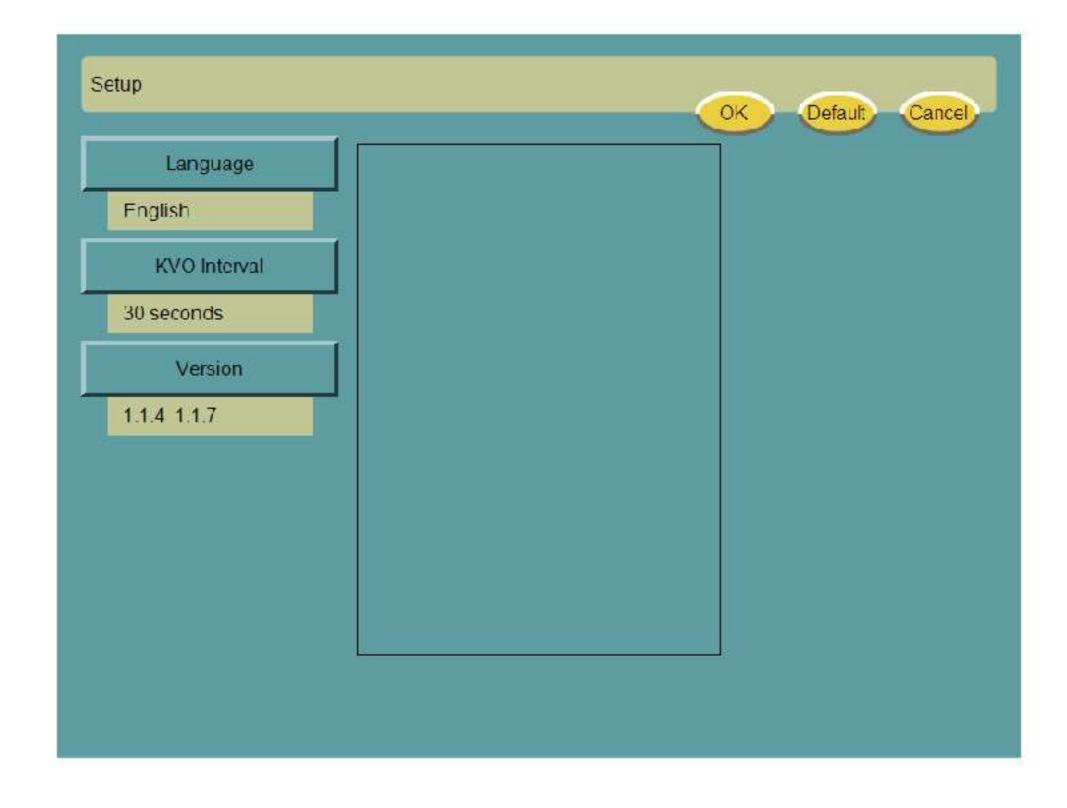

Safety

- To ensure the safety of the injection, there is a double confirmation feature with separate prompts to confirm air purging and injection ready status before any injection
- · Automatically slows down the injection when approaching to pressure limit
- · Injection stops when pressure limit is reached. Alerts with both light and sound
- Waterproof design inside the head to minimize damage from contrast media and saline leakage